

Lab 2: CDOM absorption**7 July 2015**

LABORATORY SAFETY ISSUES – isopropyl alcohol for cleaning ac-meters; general laboratory safety

INTRODUCTION

The major absorbers in seawater are **water** itself, chromophoric or color-absorbing dissolved organic matter (**CDOM**; in older literature, the term ‘g’ for Gelbstoff was used; the British used Gilvin; the German’s used yellow substance), and **absorbing particles**. This lab introduces Beer’s and its application in a standard spectrophotometer (the workhorse of chemistry) and the WET Labs absorption and attenuation meters (ac-9 and ac-s).

In the spectrophotometer lab (Station 1) the output measured is *A*, absorbance, which is unitless and based on log 10. In the in situ absorption meter lab (Station 2) the output is the absorption coefficient, *a*, with units of m^{-1} and is based on the natural log. Typically *a* is reported as the spectral absorption coefficient with the designation (λ). Please see Appendix 1 for review of Beer’s Law. Note the conversion between *A* (log 10) and *a* (natural log) is:

$$a = A \cdot (0.434 \cdot \ell)^{-1} = 2.304 \cdot A \cdot \ell^{-1}$$

where ‘ ℓ ’ is the pathlength of the instrument. The spectral absorption coefficient of CDOM ($a_{\text{CDOM}}(\lambda)$) is operationally defined as the absorption of seawater or freshwater that has been passed through a filter MINUS the absorption of a high quality water blank (such as Milli-Q type 1 water with a UV oxidizing cartridge <http://www.millipore.com/lab_water/clw4/type1>). Typically a 0.2 μm plastic filter such as a Sartoris filter or a glass fiber filter such as a Whatman G/FF filter with a nominal pore size of 0.7 μm is used.

CDOM analyses should be carried out as soon as possible after water collection and filtration because colloid formation can continue after filtration. Typically samples are collected in amber glass bottle. Sampling and processing containers should be clean.

Several things to consider in the filtration and analysis of the CDOM sample are:

- inclusion of viruses and small bacteria in the filtered fraction,
- colloidal nature of the material passing through the filter (how does a colloid form, and could colloids form while you are processing?),
- effective pore size of the filter (0.2 μm vs. 0.7 μm),
- does the effective pore size of the filter change as a function of filter loading?
- role of salts and colloidal-size particles in scattering,
- quality of the ‘pure’ water blank (how pure is pure? how do you know it is pure?),
- chemical and absorptive nature of the dissolved organic matter (DOM) and adsorbed minerals.

STATIONS:

- 1 – measure **A** (absorbance, unitless) with a spectrophotometer, Beer’s Law in a tank.
- 2 – measure **a** (absorption coefficient, units of m^{-1}) with ac-9 or ac-s.

WATER SAMPLES: water samples will be filtered through a nominally 0.2 μm water. Water samples should be at room temperature. Record the temperature.

(Save water in white buckets for fluorescence lab on Friday; store in dark).

1. DRE - Damariscotta River Estuary water (collected at dock)
2. Biscay Pond (a freshwater 'tea' lake)

Half the class goes to the Perry lab (Station 1), half to the Mitchell classroom (Station 2).

STUDENTS divide into six groups of 3 or 4 students per group.

After ~ 2 hours, students switch stations.

Station 1 – Perry Lab			
Group A	spectrophotometer with 1 cm cuvette	1) Milli-Q water	Beer's law in a tank; one group at a time
Group B	spectrophotometer with 5 cm cuvette	2) RO water	
Group C	spectrophotometer with 10 cm cuvette	3) tap water 4) 0.2 μm -filtered DRE water 5) 0.2 μm -filtered Biscay Pond water	
Station 2 – Mitchell Classroom			
Group A	ac-s	0.2 μm -filtered DRE water	
Group B	ac-9	0.2 μm -filtered DRE water	
Group C	ac-s	0.2 μm -filtered Biscay Pond water	

STATION 1 – measure A (absorbance with unitless dimensions) with a bench-top spectrophotometer (Cary-50).

Three groups; each group will perform the same measurements but will use a different pathlength cuvette (1 cm, 5 cm or 10 cm); collectively as a class, students will pool data to examine effect of pathlength (L) on absorbance (A).

1. Cary-50 spectrophotometer scan limits 200 – 800 nm & medium scan speed;
2. cuvettes have been pre-cleaned with RBS detergent:
 - don't touch optical surfaces; wipe optical surface with lens paper, NOT Kimwipes;
 - make sure all water samples are at room temperature;
 - rinse cuvette 3 times with a few mL of sample to remove previous sample;
 - visually inspect cell to ensure that there are no visual in-homogeneities (bubbles, residual mixing/turbulence between fresh and salt water sample, particles; grease);
 - place cell in holder in **EXACTLY** same orientation every time (dot on cuvette faces same direction, larger cuvettes are tipped in the same direction).

STATION 2 – measure a (absorption with units of m^{-1}) with ac-meter.

Three groups; each group will perform the same measurements, either on the ac-9 or ac-s. All instruments have 25-cm pathlengths.

Clean the sensor windows and tubes prior to measurements with lens paper and ethanol. Measure the temperature of every sample. If salinity of the Damariscotta River Estuary is not available, assume a salinity of 28.

➔ Focus on getting good pure water calibrations: each student should run her/his own Milli-Q water cal (either a-tube or c-tube of the ac-meter, or both). Save files in your group's folder.

Run the filtered water samples in both the a-tube AND c-tube of the ac-meter. Remember to save the files in your group's folder:

Group A – ac-s with 0.2 μm -filtered DRE water

Group B – ac-9 with 0.2 μm -filtered DRE water

Group C – ac-s with 0.2 μm -filtered Biscay Pond water

Instructions and code for processing data will be provided at the lab.

ASSIGNMENTS – DIVIDE THE WORK AND CONQUER! Come prepared to deliver a briefing tomorrow morning (feel free to reorganize within/among groups as you see fit).

BEER’S LAW:

The tank is an analogue spectrophotometer with light source, sample (water with green dye) and detector. Unlike the Cary 50 spectrophotometer, the light is not strictly monochromatic.

Plot I_n/I_0 vs. volume of dye.

SPECTROMETER BLANKS:

All spectrometers (Cary-50, ac-meters, etc.) depend on a pure water blank.

- 1) For the Cary-50, examine the RO and tap water spectra – would RO and tap water be sufficient as a blank? Is that true for all wavelengths?
- 2) For the ac-meters, you used Milli-Q water as the blank.
 - how does your blank compare to the device file?
 - why could your blank be different from the device file (factory blank)?
 - what are the symptoms of a bad ‘pure’ water calibration?
 - what, in addition to bad water, could contribute to a poor calibration?

TEMPERATURE EFFECT

- 1) For the ac-meters, calculate the spectral absorption coefficient, $a_{\text{CDOM}}(\lambda)$, of DRE and Biscay Pond waters for both the absorption flow tube (a tube) and the attenuation flow tube (c tube). Apply:
 - Milli-Q pure water calibration
 - corrections for temperature (at the time of measurement) and salinity (Biscay Pond is freshwater and assume DRE has a salinity of 28)
 - how would your results change if you used a temperature correction that was 2° C too high or too low? Is there a diagnostic pattern that suggests an incorrect temperature correction?
 - how would your results change if you used a salinity correction that was 1 unit too high or too low? 5 units too high or too low?
 - is there any difference between the $a_{\text{CDOM}}(\lambda)$ for the a tube vs. the c tube? Are these consistently different between or among samples?
- 2) Convert 10 cm Cary-50 absorbance values to absorption coefficients for 705 to 725 nm: $a_{\text{CDOM}}(\lambda) = 2.304 \cdot A(\lambda) \cdot \ell^{-1}$. Compare corrected ac-meter and Cary-50 data for the field samples.
 - are these values equal to zero?

- is there justification for forcing these values to zero, and subtracting the average from all other wavelengths? See below on spectral slopes (are they sensitive to NIR values?).

PATHLENGTH:

- 1) Cary-50: For each of the four groups use both DRE and Biscay Pond waters to plot absorbance (**A**) at 412 nm vs. pathlength (1, 5, 10 cm).
Is **A** a linear function of path length (Beer's Law)?
- 2) Convert ac-meter absorption coefficient at 412 nm to absorbance (natural log to log base 10): [$a_{\text{CDOM}}(\lambda) * 100 \text{ cm}$]/ 0.434 = $A(\lambda)$. Add these data to the plot of **A** vs. pathlength.
Is **A** a linear function of pathlength?

SPECTRAL SLOPE:

Calculate the spectral slope of the CDOM absorption coefficient, S_{CDOM} , for DRE and Biscay Pond water. For the ac-meters, the slope will be for the visible range only. Convert spectrophotometer values of **A** to **a**; depending on the magnitude of **A**, use the appropriate pathlength. Questions to consider are:

- using data from the Cary-50, is the spectral slope constant as a function of wavelength?
 - what is the slope for the UV? Does it change for different regions within the UV?
 - what is the slope in the visible for the range 412 nm – 676 nm for the ac-meters and for the spectrophotometer?
- compare spectral slopes calculated using data from different instruments.
- does the spectral slope vary as a function of the method used to compute the slope?
- does the slope change if wavelengths above 700 nm are omitted, included?
- does the spectral slope vary as a function of water source?

Equation for and calculation of spectral slope:

$$a_{\text{CDOM}}(\lambda) = a_{\text{CDOM}}(\lambda_{\text{REF}}) e^{-S_{\text{CDOM}}(\lambda - \lambda_{\text{REF}})}$$

where $a_{\text{CDOM}}(\lambda)$ is the amplitude of the absorption coefficient at any wavelength λ (Jerlov, 1976) or at the reference wavelength, λ_{REF} (usually 412 or 440 nm). See Carder et al. (1989) and Blough and Del Vecchio (2002) for a discussion of the interpretation of the spectral slope.

The best method to calculate the slope is to minimize the square difference between the exponential model and the data (possibly weighed by a different error in each


wavelength if the uncertainty varies as function of wavelength, e.g. due to variability in source intensity as function of wavelength). The relative (percent) error is not constant spectrally; in the red the absorption is low and the signal-to-noise high. Slope measurements often exclude red wavelengths due to its sensitivity to temperature (e.g. the 715 nm channel in the ac-9).

You may write your own code to determine the slope by non-linear exponential regression (we will also supply code: <http://misclab.umeoce.maine.edu/software.php>); OR,

less rigorously, you may determine the spectral slope for $a_{\text{CDOM}}(\lambda)$ by plotting the ln-transformed values of a_{CDOM} vs. wavelength using Excel and adding a trend line (this is the same as if you fit an exponential curve in Excel, try it). If you use the latter method, *is the slope linear?*

APPENDIX 1 REVIEW OF BEER'S LAW

<http://teaching.shu.ac.uk/hwb/chemistry/tutorials/molspec/beers1.htm>

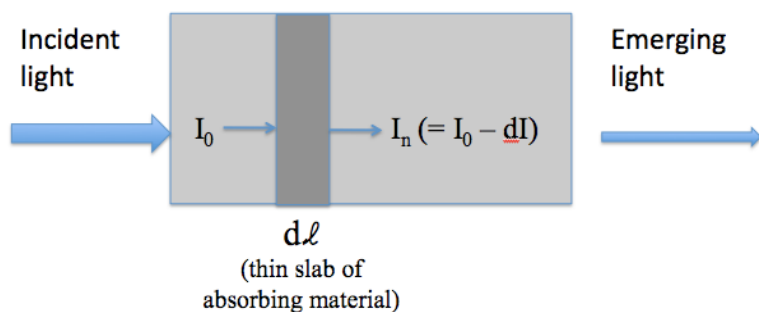

$$-dI = I_0 \cdot (\epsilon \cdot C) \cdot d\ell \quad (-dI \text{ is fraction of light absorbed in thin slab})$$

$$-dI / I_0 = \epsilon \cdot C \cdot d\ell$$

$$I_n = I_0 \exp(-\epsilon \cdot C \cdot \ell) \quad (\text{integrating over entire pathlength})$$

$$I_n / I_0 = \exp(-\epsilon \cdot C \cdot \ell) \quad (\text{also termed transmittance})$$

$$I_0 / I_n = \exp(\epsilon \cdot C \cdot \ell)$$

$$\ln I_0 / I_n = \epsilon \cdot C \cdot \ell = a \cdot \ell$$

where:

I_0 is intensity of light before it passes through the sample,

I_n is the intensity measured at the detector after light passes through the sample,

I_n / I_0 is unitless,

ϵ is the molar absorption coefficient – a measure of how much light a 1 M solution of dye will absorb ($\text{m}^2 \text{mol}^{-1}$),

C is the concentration of the dye (mol m^{-3}), and

ℓ is the path length that the light must travel through the solution (m).

NB: Here we combine terms “ $\epsilon \cdot C$ ” into a single term “ a ”, the absorption coefficient (m^{-1}).

Notice that the Beer's Law equation is written in log base e (natural logarithms, ln).

However, spectroscopists historically used log base 10, rather than log base e. The principle is the same but A , the absorbance output from the spectrophotometer, is log base 10. Also, chemists include pathlength in A , because all measurements are typically made using same pathlength; hence, A is reported with as dimensionless and the value of A will change with pathlength.

$$\begin{aligned} A &= \log_{10}(I_0/I_n) \\ &= \log_{10}e \cdot \ln(I_0/I_n) \\ &= 0.434 \cdot \ln(I_0/I_n) = 0.434 \cdot (a \cdot \ell) \end{aligned}$$

$$a = A \cdot (0.434 \cdot \ell)^{-1} = 2.304 \cdot A \cdot \ell^{-1}$$

Remember from calculus, that when changing log bases: $\log_a X = \log_b X \cdot \log_a b$.

To convert a natural logarithm to a base 10 logarithm, multiply by $\log_{10}e (=0.434)$.

To convert a base 10 logarithm to a natural logarithm, multiply by $\log_e 10 (=2.304)$.

Limitations of the Beer-Lambert law

The linearity of the Beer-Lambert law is limited by chemical and instrumental factors.

Causes of nonlinearity include:

- deviations in absorptivity coefficients at high concentrations due to electrostatic interactions between molecules in close proximity
- pathlength amplification due to scattering of light by particulates in the sample
- fluorescence or phosphorescence of the sample
- changes in refractive index at high analyte concentration
- shifts in chemical equilibria as a function of concentration
- non-monochromatic radiation, deviations can be minimized by using a relatively flat part of the absorption spectrum such as the maximum of an absorption band
- stray light

Source: <http://www.chemistry.adelaide.edu.au/external/soc-rel/content/beerslaw.htm>

BACKGROUND MATERIAL ON SPECTROPHOTOMETRY

General principles of operation of spectrophotometer:

Web sites that present good reference material on the fundamentals of UV-visible spectrometry:

<http://www.cem.msu.edu/~reusch/VirtualText/Spectrpy/UV-Vis/uvspec.htm#uv1>

<http://www.cem.msu.edu/~reusch/VirtualText/Spectrpy/InfraRed/infrared.htm>

Across the Spectrum: Instrumentation for UV/Vis Spectrophotometry

Slightly modified and shortened from Shane Beck, 1998, *The Scientist*, 12(3): 20.

Modern spectrophotometry was pioneered by Dr. Arnold Beckman in the 1940's.


1. Light source: typical UV/Vis spectrophotometers utilize two light sources: a deuterium arc lamp for consistent intensity in the UV range (190 to 380 nm) and a tungsten-halogen lamp for consistent intensity in the visible spectrum (380 to about 800 nm). Some spectrophotometers, such as the Cary 50, have a xenon flash lamp.

2. Dispersion of light into different wavelengths can occur before or after the light passes through the sample. The monochromator disperses light into different angles by prisms or holographic gratings. NB: with a prism, the angle of dispersion can be nonlinear and sensitive to changes in temperature. In contrast, holographic gratings eliminate nonlinear dispersion and are not temperature sensitive; they are glass blanks with narrow ruled grooves. The grating itself is usually coated with aluminum to create a reflecting source. Gratings do require filters since light is reflected in different orders with overlapping wavelengths.


Light passing through the monochromator exits as a band. The width of this band of light at half the maximum intensity is the spectral bandwidth. Bandwidth comes in to play with regard to accuracy, since the accuracy of any absorbance measurement is dependent on the ratio of the spectral bandwidth to the natural bandwidth of the substance being measured. The natural bandwidth is the width of the absorption band of the sample at half the absorption maximum. As a rule, a ratio between spectral bandwidth and natural

bandwidth of 0.1 or less will generate absorbance measurements 99.5 percent accurate or better. Above this, accuracy deteriorates.

Adapted from Hewlett-Packard's "Fundamentals of Modern UV-Visible Spectroscopy" publication number 12-5965-512E, 1996.

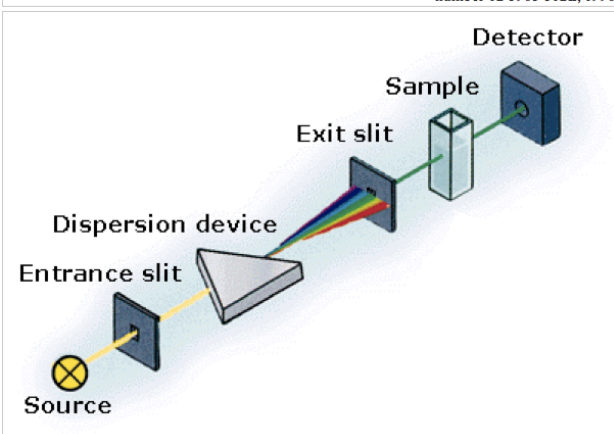

Adapted from Hewlett-Packard's "Fundamentals of Modern UV-Visible Spectroscopy" publication number 12-5965-512E, 1996.

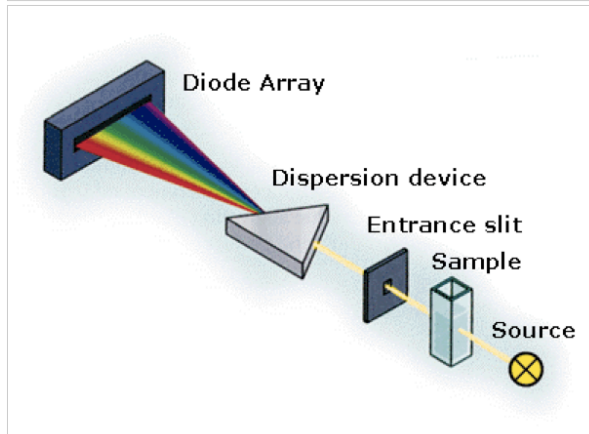

4. Sample absorbance is determined by comparing the intensity of the light passing through the sample and hitting the detector vs. intensity of light passing through a blank. Detectors include: a) photomultiplier tube, with good sensitivity throughout the UV/Visible spectral range and highly sensitive at low light levels; or b) photodiode with a wider dynamic range, and consisting of a semiconductor and a capacitor to charge the semiconductor. As light hits the semiconductor, electrons flow through it, thereby lowering the charge on the capacitor. The intensity of light of the sample is proportional to the amount of charge needed to recharge the capacitor at predetermined intervals. Often the detector is composed of a photodiode array, with photodiode detectors arranged on a silicon crystal so a spectral scan is instantaneous.

In single-beam spectrophotometers, the blank and sample are not measured simultaneously. Interspersing measurements of samples and blanks are needed to correct for lamp drift. Dual-beam spectrophotometers utilize a "chopper" or beam splitter that alternates the light path between the reference optical path and sample optical path to the detector at a speed that minimizes medium- or long-term effects of lamp drift. Some dual beam instruments scan continuously so that the sample, blank and dark reference are actually performed at different wavelengths (leading to a skewing effect as a function of wavelength, dependent upon scan speed); in others there is a phase locked wavelength drive so that the sample, blank and dark reference readings occur at the same wavelength. Blanks should be refreshed to prevent sample warming, or kept in a cooled holder.

If the sample is not a pure solution, scattering can occur. An integrating sphere can be used to collect scattered light, and correct the instrument reading to provide true absorption. The coatings on integrating spheres are highly scattering so as to ensure that the light field within the sphere is isotropic and therefore measuring a small portion of that light is equivalent to measuring it all. However, the coatings are also particularly absorptive of UV and blue radiation, which limits their utility in the UV range.

APPENDIX 1 - BACKGROUND MATERIAL ON IN SITU SPECTROPHOTOMETERS

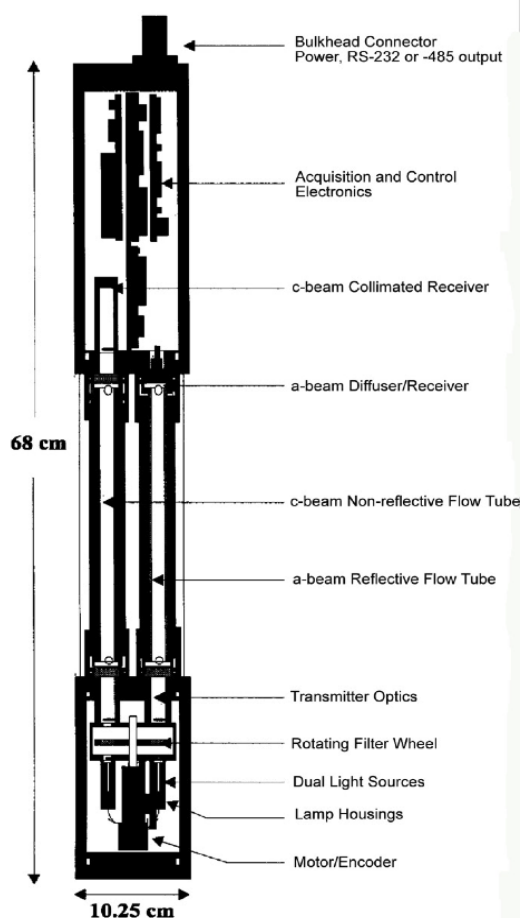


Figure 1. General schematics of WET Labs ac9.

General principles of operation of the ac-9:

The only commercially available mature *in situ* absorption meter is manufactured by WET Labs (Figure 1) <http://www.wetlabs.com>.

Some important issues related to using the ac-9 can be found in Pegau et al. 1995; Pegau et al. 1997; Bricaud et al. 1995; Zaneveld et al (1994); Twardowski et al., 1999; Roesler and Boss (2007); Leymarie et al., (2010).

Schematics (Figure 2)

1. Light source: Incandescent bulb
2. Dispersion of light into different wavelengths is done by a filter wheel with 9 filters. The filter wheel spins at 6 Hz yielding 6 spectra per second. The filter band width is 10 nm.
3. A collimated beam of light passes through the sample and onto a diffuser and a single diode detector (in the case of *a*; to maximize the capture of forward scattered light) or into a narrow angle detector (in the case of *c*, to minimize the capture of forward scattered light).
4. Sample absorption is determined relative to a pure water calibration provided by the factory (contained in the device file, ac-90nnn.dev, where nnn is the instrument serial number). Given the tendency for drift and alignment issues, it is standard practice to run your own pure water calibration prior and subsequent to each experiment.

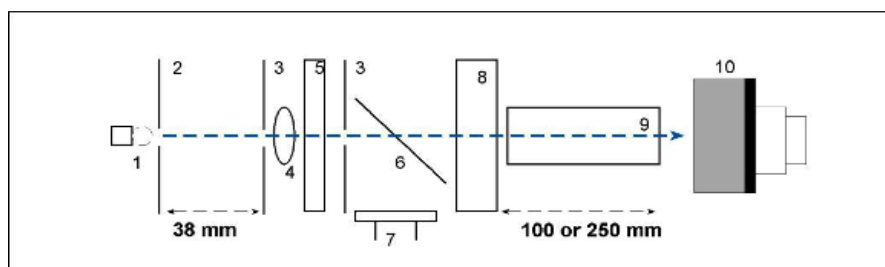


Figure 2. Schematic Representation of absorption beam optics

- 1 Lamp
- 2 1 mm aperture
- 3 6 mm aperture
- 4 38 mm singlet lens
- 5 Interference filter

- 6 Beam splitter
- 7 Reference detector
- 8 6 mm quartz pressure window
- 9 Reflective flow tube
- 10 Diffuser/Signal detector



Figure 3. Filter wheel of acs.

The hyperspectral version of the ac-9 is called the ac-s. Although similar in design the filter wheel holds two sections of a Linear Variable Filter (LVF), centered 180 degrees from each other on the filter wheel (Figure 3). The two filter sections are cut from a single LVF such that a portion of the spectrum around 550 nm is covered by both filters. This overlap is to allow for merging of the data from both filter sections (data generally display a slight error at this overlap that needs to be corrected for). Each filter covers approximately a 72 degree

section of the beam path across the filter wheel. The filter wheel rotates at a tightly controlled 8.0 rps, such that the shorter wavelength of each filter section is traversed before the longer wavelength.

REFERENCES

- Babin, M., D. Stramski, G.M. Ferrari, H. Claustre, A. Bricaud, G. Obolensky, and N. Hoepffner. 2003. Variations in the light absorption coefficients of phytoplankton, nonalgal particles, and dissolved organic matter in coastal waters around Europe. *Journal of Geophysical Research Ocean* 108 (C7): article number 3211 , doi: 10.1029/2001JC000882.
- Blough, N.V., and R. Del Vecchio. 2002. Chromophoric dissolved organic matter (CDOM) in the coastal environment. In: D. Hansell and C. Carlson, Editors, *Biogeochemistry of Marine Dissolved Organic Matter*, Academic Press, San Diego, CA.
- Bricaud, A., A. Morel, and L. Prieur. 1981. Absorption by dissolved organic matter of the sea (yellow substance) in the UV and visible domains. *Limnol. Oceanogr.* 26: 43-53.
- Carder, K. L., R. G. Steward, G. R. Harvey, and P. B. Ortner. 1989. Marine humic and fulvic acids: Their effects on remote sensing of ocean chlorophyll. *Limnol. Oceanogr.* 34: 68-81.
- Jerlov, N.G. 1976. *Marine Optics*. Elsevier, New York.
- Leymarie, E., D. Doxaran, and M. Babin. 2010. Uncertainties associated to measurements of inherent optical properties in natural waters. *Appl. Opt.* 49(28): 5415-5436.
- Pegau W. S., G. Deric and J. R. V. Zaneveld. 1997. Absorption and attenuation of visible and near-infrared light in water: dependence on temperature and salinity. *Applied Optics.* 36: 6035-6046.
- Pegau, W. S., J. S. Cleveland, W. Doss, C. D. Kennedy, R. A. Maffione, J. L. Mueller, R. Stone, C. C. Trees, A. D. Weidemann, W. H. Wells and J.R.V. Zaneveld. 1995. A comparison of methods for the measurement of the absorption coefficient in natural. *J. Geophys. Res.* 100(C7): 13201-13220.

- Roesler, C.S., 1998: Theoretical and experimental approaches to improve the accuracy of particulate absorption coefficients derived from the quantitative filter technique. *Limnology and Oceanography*. 43: 1,649-1,660.
- Roesler, C. S. and E. Boss. 2007. In situ measurement of the inherent optical properties (IOPs) and potential for harmful algal bloom (HAB) detection and coastal ecosystem observations, Chapter 5, pp. 153-206. In Babin, M, C. S. Roesler and J. J. Cullen [eds.] *Real-Time Coastal Observing Systems for Ecosystem Dynamics and Harmful Algal Blooms*. UNESCO Series Monographs on Oceanographic Methodology.
- Twardowski M. S., J. M. Sullivan, P. L. Donaghay and J. R. V. Zaneveld. 1999. Microscale quantification of the absorption by dissolved and particulate material in coastal waters with an ac-9. *Journal of Atmospheric and Oceanic Technology*. 16: 691-707.
- Zaneveld J. R. V., J. C. Kitchen and C. M. Moore. 1994. The scattering error correction of reflecting-tube absorption meters. *Proc. SPIE, Ocean Optics XII*. 2258: 44-55.