



***ECO* Fluorometers and Scattering Sensors**

User manual

05/2014, Edition 3



Section 1 Specifications	3
1.1 Mechanical.....	3
1.1.1 Bulkhead connectors.....	3
1.2 Electrical.....	4
1.3 Communications.....	4
1.4 Optical.....	4
1.4.1 Single-parameter fluorometer.....	4
1.4.2 Single-parameter scattering.....	4
1.4.3 Two-parameter fluorometer-turbidity.....	4
1.4.4 Three-parameter fluorometer/scattering.....	5
Section 2 Operation and maintenance	7
2.1 Check sensor operation.....	7
2.2 Set up sensor for deployment.....	8
2.3 Monitor data.....	8
2.3.1 Monitor engineering units output.....	8
2.4 Get data from sensor.....	9
2.5 Additional operations.....	10
2.5.1 Set date and time.....	10
2.5.2 Adjust option to collect data.....	10
2.5.3 Adjust views in plot data tab.....	11
2.6 Sensor maintenance.....	12
2.6.1 Bulkhead connector maintenance.....	12
Section 3 Reference topics	13
3.1 Delivered items.....	13
3.2 Calibration.....	13
3.3 Characterization.....	13
3.4 Field characterization.....	13
3.4.1 Store field characterization values in device file.....	14
3.4.2 Store field characterization values in sensor.....	14
3.5 Chlorophyll fluorometer-specific setup.....	15
3.6 Device files.....	15
3.7 Terminal program operation.....	17
3.7.1 Common terminal program operations.....	17
3.7.2 Terminal program operation for sensors with internal memory.....	17
3.7.3 ECOView and terminal program differences.....	18
Section 4 Optional equipment	19
4.1 Test cable.....	19
4.2 Internal batteries.....	19
4.2.1 Replace batteries.....	19
4.3 Bio-wiper and copper faceplate.....	22
4.3.1 Clean Bio-wiper and faceplate.....	22
4.4 External thermistor.....	24
4.5 Pressure sensor.....	24
4.5.1 Maintain pressure sensor.....	25
Section 5 General information	27
5.1 Warranty.....	27
5.2 Service and repair.....	27
5.3 Support email addresses.....	27
5.4 Waste electrical and electronic equipment.....	27
Appendix A ECO mounting bracket	29

Section 1 Specifications

ECO sensors measure different parameters of the Earth's natural waters. The sensors are available in several models and with a variety of optional features.


Real-time (RT)	Output is digital or analog. Operates continuously. Does not store data.
Real-time deep (RTD)	6,000 m rated depth. Does not store data.
Standard	Output is digital or analog. Has low power mode. Stores data.
Bio-wiper (S)	Standard, and with Bio-wiper to prevent bio-fouling.
Battery (B)	Standard, and with internal batteries.
Bio-wiper and battery (SB)	Standard, with both a Bio-wiper and internal batteries.


Note: The three-parameter sensor is available in Standard and Battery models only and has no analog output.

1.1 Mechanical

	RT, standard	RTD	S	B, SB
Diameter	6.30 cm			
Length	12.70 cm	17.68 cm	13.3 cm	28.00 cm
Depth rating	600 m	6000 m	300 m	
Temperature range	0–30 °C			
Weight in air, water	0.40 kg, 0.02 kg	1.3 kg, 0.75 kg	0.50 kg, 0.08 kg	0.96 kg, 0.14 kg

1.1.1 Bulkhead connectors

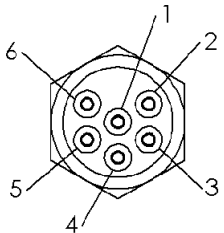Pin	Function	MCBH-6-MP
1	Ground	
2	RX	
3	Analog 2	
4	Voltage in	
5	TX	
6	Analog out	

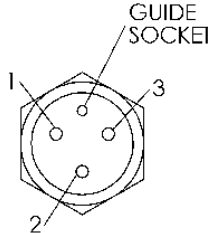Additional bulkhead connector on sensors with internal batteries.		
Socket	Function	Diagram of MCBH-3-FS
1	Voltage in	
2	No connect	
3	Battery out	

Note that the three-pin blue-tipped jumper plug supplied with the sensor can be used to supply power to the sensor.

Specifications

1.2 Electrical

Input	7–15 VDC
Current draw, typical	50 mA; 60 mA (triplet)
Current draw, standby	140 μ A
Current draw, active wiper	140 mA
Linearity	99%

1.3 Communications

	Fluorometer	NTU	Scattering	FLNTU	Triples
Sample rate	to 8 Hz				to 4 Hz
Data storage	108000 samples			90000 samples	77000 samples
RS232 output rate	19200 baud				
Data resolution	14 bit		12 bit		
Digital output maximum	~16380 counts		4130 \pm 30 counts		
Analog output maximum	5 V				no analog

1.4 Optical

1.4.1 Single-parameter fluorometer

Parameter	Wavelength EX/EM	Range, Sensitivity
Chlorophyll (Chl)	470/695 nm	0–125, 0.016 μ g/L
Colored Dissolved Organic Matter (CDOM)	370/460 nm	0–500, 0.093 ppb
Uranine (UR)	470/530 nm	0–400, 0.05 ppb
Phycocyanin (PC)	630/680 nm	0–230, 0.029 ppb
Phycoerythrin (PE)	540/570 nm	0–230, 0.029 ppb

1.4.2 Single-parameter scattering

Parameter	Wavelength	Range, Sensitivity
Scattering	470 nm, 532 nm, 650 nm	0–5, 0.003 m^{-1}
	700 nm	0–3, 0.002 m^{-1}
		0–5, 0.003 m^{-1}

1.4.3 Two-parameter fluorometer-turbidity

Parameter	Wavelength EX/EM	Range, Sensitivity (chl)	Parameter	Wavelength	Range, Sensitivity (NTU)
Chlorophyll	470/695 nm	0–30, 0.015 μ g/L 0–50, 0.025 μ g/L 0–75, 0.037 μ g/L 0–125, 0.062 μ g/L 0–250, 0.123 μ g/L	NTU	700 nm	0–10, 0.005 NTU 0–25, 0.013 NTU 0–200, 0.098 NTU 0–350, 0.172 NTU 0–1000, 0.123 NTU

1.4.4 Three-parameter fluorometer/scattering

Parameter	Wavelength EX/EM	Range, Sensitivity
Chlorophyll (Chl)	470/695 nm	0–30, 0.015 µg/L
		0–50, 0.025 µg/L
Colored Dissolved Organic Matter (CDOM)	370/460 nm	0–375, 0.184 ppb
Uranine (UR)	470/530 nm	0–300, 0.073 ppb
Phycocyanin (PC)	630/680 nm	0–175, 0.086 ppb
Phycoerythrin (PE)	540/570 nm	0–175, 0.086 ppb

Parameter	Wavelength	Range, Sensitivity
Scattering	412 nm, 470 nm, 532 nm, 650 nm, 880 nm	0–5, 0.003 m ⁻¹
	700 nm	0–3, 0.002 m ⁻¹
		0–5, 0.003 m ⁻¹

Section 2 Operation and maintenance


2.1 Check sensor operation

⚠ CAUTION

CDOM sensors use UV LED light. Do not look directly at a UV LED when it is on. It can damage the eyes. Keep products that have UV LEDs away from children, pets, and other living organisms. Wear polycarbonate UV-resistant safety glasses to protect the eyes when a UV LED is on.

Make sure that the sensor operates before further setup and deployment.

1. Connect the 6-socket connector on the optional test cable (refer to the section on the [Test cable](#) on page 19 for details) to the sensor.
2. Remove the cap that protects the sensor's optical face.
3. Connect a serial-to-USB adapter to the test cable to connect the test cable to the host PC.
4. Connect the sensor to a power supply.
 - a. Connect sensors with internal batteries to the manufacturer-supplied, three-pin, blue-tipped jumper plug.
The sensor comes on.
 - b. Connect sensors without internal batteries to the optional test cable and a regulated power supply set at 12 VDC.
5. Start the host software from the manufacturer-supplied CD.
 - a. Select the COM port on the host PC.
 - b. Select the device file for the sensor from the CD.
 - c. Select the baud rate if necessary. The default is 19200.

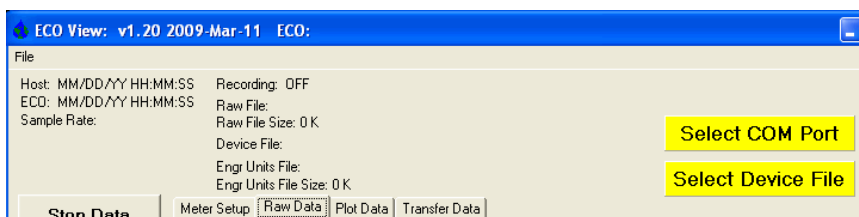

6. Turn on the power supply.
The sensor comes on.
7. Push **Start Data** on the host software.
8. Go to the *Raw Data* tab of the host software.
The data collected by the sensor shows in the "Signal" column.

Figure 1 Format of data collected by most ECO sensors

Meter Setup		Raw Data	Plot Data	Transfer Data			
06/14/12 06:06:05	532	267	660	3070	695	78	535
06/14/12 06:06:06	532	315	660	3406	695	89	535
06/14/12 06:06:07	532	437	660	3861	695	127	535
06/14/12 06:06:08	532	509	660	4122	695	175	535
06/14/12 06:06:09	532	2577	660	4122	695	629	535
06/14/12 06:06:10	532	4122	660	4122	695	957	535
06/14/12 06:06:11	532	4122	660	4122	695	970	535
06/14/12 06:06:12	532	4122	660	4122	695	868	535
06/14/12 06:06:14	532	4122	660	4122	695	946	535

Date Time Wave-length Signal Wave-length Signal Wave-length Signal Thermistor

Note that RT and puck sensors show 9's instead of the Date and the Time.

Figure 2 Format of data collected by real-time sensors

99/99/99	99:99:99	695	42	700	264	460	51	538
99/99/99	99:99:99	695	43	700	260	460	55	538
99/99/99	99:99:99	695	41	700	257	460	64	538
99/99/99	99:99:99	695	37	700	255	460	62	538
99/99/99	99:99:99	695	39	700	258	460	50	538
99/99/99	99:99:99	695	44	700	262	460	53	538
99/99/99	99:99:99	695	49	700	259	460	58	538
	Wave-length	Signal	Wave-length	Signal	Wave-length	Signal	Thermistor	

- Look at the maximum data value for the sensor. Hold a finger, the protective cap, or fluorescent stick if the sensor is a fluorometer, 1–4 cm away from the sensor's optics face.
The data value in the "Signal" column in the *Raw Data* tab will increase toward the maximum data value specified for the sensor.
 - Scattering and turbidity sensors: use a finger or the protective cap.
 - CDOM fluorescence sensors: use the blue fluorescent stick.
 - Chlorophyll or phycoerythrin sensors: use the orange fluorescent stick.
 - Uranine or phycocyanin sensors: use the yellow fluorescent stick.
 - PAR sensors: point the sensor into the light.
- Push **Stop Data**.
The Bio-wiper closes on sensors so-equipped. If the power is turned off in mid-cycle, the Bio-wiper starts at the beginning of the cycle when power is supplied again.

2.2 Set up sensor for deployment

- Refer to the previous section to make sure that the sensor operates correctly.
- Replace the test cable with a sea cable for deployment.
- Remove the protective cap from the sensor if necessary.
- Use either the blue-tipped jumper plug or an external power supply to supply power to the sensor for deployment.
- Refer to the section on [Additional operations](#) on page 10 for details on how to set up the sensor for a given application.

2.3 Monitor data


Monitor the data from the sensor in counts. The number of "Signal" columns will vary depending on whether the user has a one-, two-, or three-parameter sensor.

- Make sure the sensor has power supplied and is turned on.
- Push **Start Data**.
- Go to the *Raw Data* tab.
Refer to the [Check sensor operation](#) on page 7 illustration to see the format of the collected data.


Note: *RT and puck sensors usually show 9's as placeholders in the date and time columns.*

2.3.1 Monitor engineering units output

- Go to the *Plot Data* tab.
- Select "Engr Units" from the drop-down menu near the top of the tab.

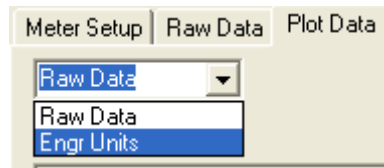

3. Select the type of units to see.

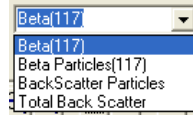

The host software calculates the engineering units for the *Plot Data* tab.

Note: Data is saved in counts, not engineering units.


2.4 Get data from sensor

⚠ CAUTION


Do data transfers away from harsh environments such as strong electric fields or electrostatic discharge sources. Electrostatic Discharge (ESD) sources may temporarily disrupt data transfer. If this occurs, move the sensor away from the ESD source. Turn the power off and then on and continue operation.

Get the data from the sensor to the host PC.

1. Make sure that the sensor is not on but has power supplied to it.
2. Select the *Transfer Data* tab in the host software.
3. Push **Receive Data**.

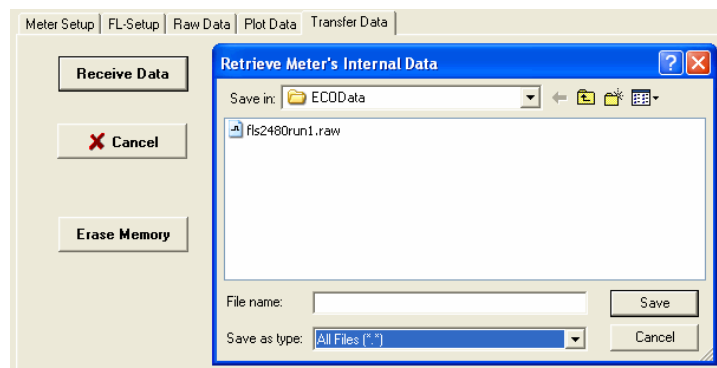

4. Select a location on the host PC to store the data.
5. Type a file name in the *Retrieve Meter's Internal Data* window.
6. Push **Save**.
The host software saves the sensor's internal data to the host PC.
7. Make sure that the data move is complete.

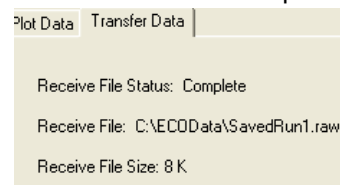

8. Open the data file to make sure the data is on the host PC.
Contact the manufacturer to get a spreadsheet template for an *ECO* sensor.
9. To erase the data from the sensor's memory, Push **Erase Memory**.

2.5 Additional operations

Note: Real-time (RT and RTD) sensors do not store data. Some of the data collection options in this section do not apply to these sensor models.


2.5.1 Set date and time

Make sure that the sensor is connected to a power supply that is turned on. Make sure that the host software program is open.

1. If the sensor is in operation, push **Stop Data** to stop the sensor.
2. Push **Set Date and Time** in the host software.
The host software adjusts the sensor's time to agree with the host PC.



3. Push **Get Date/Time/Setup** to make sure that the sensor and the host PC show the same current time.

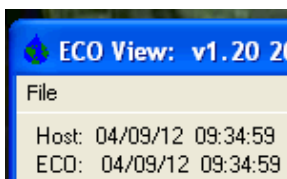


2.5.2 Adjust option to collect data

The manufacturer sets *ECO* sensors to operate at approximately 1 Hz, with the optional internal data storage turned on.

Table 1 Options to collect data

ECOView option	How it works
Set Avg/Data Rate	Adjust from 1 to 65535. Examples:
	1-parameter sensors—approximately 1 Hz = 65; approximately 2 Hz = 30
	2-parameter sensors—approximately 1 Hz = 30; approximately 2 Hz = 15
	3-parameter sensors—approximately 1 Hz = 18; approximately 2 Hz = 6
	PAR sensors—approximately 1 Hz = 310; approximately 2 Hz = 170
Set Number of Samples	Adjust from 0 to 65535. Use 0 to operate continuously.
The three options below apply only to sensors equipped with internal memory	
Set Number of Cycles	Adjust from 0 to 65535. Select the number of sample groups that the sensor will collect between low-power states.
Set Cycle Interval	Adjust the time interval between sample cycles. Do not type the colons. The minimum is 5 seconds.
Turn Logging ON	Push to turn data storage off or on (sensors with internal data storage only).

Adjust any of the options to collect data in the *Meter Setup* tab.

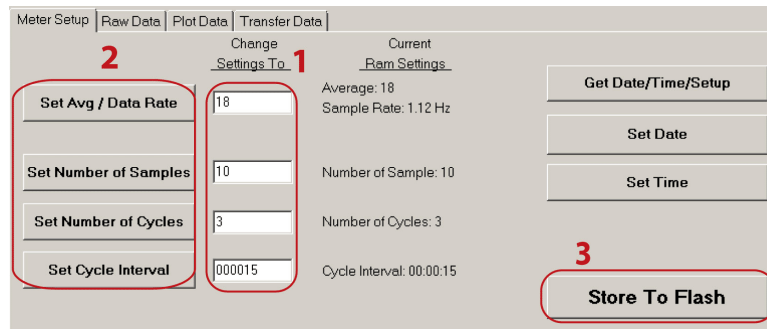

<p>1 Enter the new value into the variable box.</p>	<p>2 Push Set Avg/Data Rate that is to the left of the variable box.</p>	<p>3 Push Store to Flash. The new value will be in the column <u>Current Ram Settings</u>.</p>
--	--	--

In the example above, the 3-parameter sensor will operate at an "Average" of 18 and a "Data Rate" of 1.12 Hz. The sensor will collect data 10 rows of data for 3 cycles, with a low-power interval of 15 seconds after each cycle. The sensor stops after the last 10 rows of data are collected.

Table 2 Data collection examples

<p>Collect moored data Set Avg/Data Rate = ± 1 Hz Set Number of Samples = 50 Set Number of Cycles = 24 Set Cycle Interval = 006000 Turn Logging ON/OFF = ON The sensor will collect data once per second, 50 times every 60 minutes for 24 hours, and store the collected data.</p>	<p>Collect profiling data Set Avg/Data Rate = ± 1 Hz Set Number of Samples = 0 Set Number of Cycles = N/A Set Cycle Interval = N/A Turn Logging ON/OFF = ON The sensor will collect data once per second and store collected data until the power is turned off.</p>
---	--

If the sensor is set up to collect data intermittently, as for a moored deployment, it may be in a low-power condition. It is not possible to communicate with the sensor in this condition.

1. To start communication again, disconnect the sensor's power supply for one minute.
2. Connect the power again and push **Stop Data** several times.
3. Select the *Meter Setup* tab. Refer to [Adjust option to collect data](#) on page 10.
4. Type **0** into the **Number of Samples** variable box.
5. Push **Set Number of Samples**.
6. Push **Store to Flash**.
The sensor will now run continuously.

Make sure that the sensor operates continuously.

1. Select the *Raw Data* tab.
2. Push **Start Data**.
3. Let the sensor operate for 10 samples or more.
4. Push **Stop Data**.

2.5.3 Adjust views in plot data tab

The *Plot Data* tab in the host software lets the user see the data that the sensor has collected.



Operation and maintenance

Button	Function	Description
1	Resume	Push to start or stop the scrolling data
2	Pause	The x-axis will stop scrolling.
3	Axes scroll	Move either axis up or down, or right or left.
4	Axes zoom	Move either axis up or down, or right or left.
5	Zoom out	Decrease the detail by 2x.
6	Zoom in	Increase the detail by 2x.
7	Zoom box	Draw a box around the data area to zoom all axes.
8	Cursor	Move the cursor bar to a certain data point.
9	Copy	Copy the current data plot to the host PC clipboard.
10	Save	Save an image of the current data plot to the host PC.
11	Print	Send a snapshot of the current data plot to a printer.

View the data in either counts or in processed units ($\mu\text{g/L}$, ppb, scattering, etc.). Select the type of output from the drop-down menu above the view area.

2.6 Sensor maintenance

NOTICE

Do not use acetone or other solvents to clean the sensor or the sensor's optical face.

1. After each cast or exposure to natural water, flush the sensor with clean fresh water.
2. Use soapy water to clean any grease or oil on the sensor's optical face. It is made of ABS plastic and optical epoxy and can be damaged if an abrasive cleaner is used.
3. Dry the sensor with a clean soft cloth.

2.6.1 Bulkhead connector maintenance

Lubricate the mating surfaces of bulkhead connectors at regular intervals with pure silicone spray only. Allow the contacts to dry before they are connected.

Make sure that the pins have no corrosion, which looks green and dull. Make sure that the rubber seals on the pins are not delaminated. Connectors should mate smoothly and not feel "gritty" or too resistant.

The manufacturer recommends 3M™ Silicone Lubricant spray (UPC 021200-85822). Other silicone sprays may contain hydrocarbon solvents that damage rubber.

DO NOT use silicone grease. **DO NOT** use WD-40®. The wrong lubricant will cause the bulkhead connector to fail prematurely and the sensor will flood.

3.1 Delivered items

- the *ECO* sensor
- a dummy plug with a lock collar
- a plastic protective cover for the optics face
- a model-specific spare parts kit
- a stainless steel mounting bracket and hardware (Sensors that are depth-rated to 6000m and sensors with internal batteries do not have this bracket.)
- **On the CD:**
- this user manual
- the ECOView host software
- the sensor's device file or files
- the sensor's characterization or calibration page

3.2 Calibration

The manufacturer calibrates all scattering sensors to make sure that the data that is collected meets the sensor's specifications. This information is on the sensor-specific calibration page that comes with the sensor.

3.3 Characterization

The manufacturer characterizes all fluorescence sensors using a fluorescing material to make sure that the data that is collected meets the sensor's specifications. This information is on the sensor-specific characterization page that ships with the sensor.

3.4 Field characterization

The manufacturer recommends that the user perform a field characterization on fluorometers to make sure that the data is as accurate as possible for the user's application. The scale factor and the dark counts values can vary depending on the natural water, temperature, length of the cable, the power supply, and other factors. Do the steps below to field-characterize the sensor.

- **x** = a solution of a known concentration in volts or counts.
 - **output** = the measured sample of interest in volts or counts.
 - **dark counts** = the measured signal output in volts or counts of the sensor in clean water with black tape over the detector.
 - **scale factor** = the multiplier in $\mu\text{g/L/volt}$, ppb/L/volt , OR $\mu\text{g/L/count}$, ppb/L/count .
1. Get a solution of a known concentration, **x**.
 2. Measure and record this solution using the sensor.
This value is the **output** in volts or in counts.
 3. Measure and record the sensor's **dark counts**.
 4. Use this equation to determine the sensor's **scale factor**:
 $\text{Scale factor} = x \div (\text{output} - \text{dark counts})$.
 5. Use the scale factor to determine the concentration of the sample of interest:
 $(\text{output counts} - \text{dark counts}) \times \text{scale factor} = \text{concentration of solution}$.
 6. Store the scale factor and dark counts (offset) in the device file for the sensor, the internal memory of the sensor, or both.

3.4.1 Store field characterization values in device file

The host software uses a device file to process data. Refer to the example below.

Note: A colon comes before the comments in the device file. The comments are not used by the host software.

```

ECO          FLS-1822
Created on:  04/29/2011

:chl=ug/l.
:iengrunits=ug/l for chl; ppb for PC, PE CDOM and uranine.
:column 4=input the scale factor and offset in this column.
:N/U=not used.

maxvoltage=4.96
asv1=6.2606
asv2=12.5355
asv4=25.2860


COLUMNS=5
N/U=1
N/U=2
N/U=3
Chl=4    0.0052  48
N/U=5
    
```

1. Replace the values in column 4 of this device file with the scale factor and offset values from the field characterization.
2. Save this device file with a new name.
3. To use this new file in the host software, select the *File* menu, then push **Load Device File**.

3.4.2 Store field characterization values in sensor

Make sure that the sensor is connected to a 9V battery or power supply and the host PC to do these steps.

1. Start the host software if necessary.
2. Push **Select COM Port**. Select the communications port on the PC.
3. Push **Device File**. Select the device file from the field characterization.
4. Make sure the sensor is not collecting data.
5. Select the *FL Setup* tab.

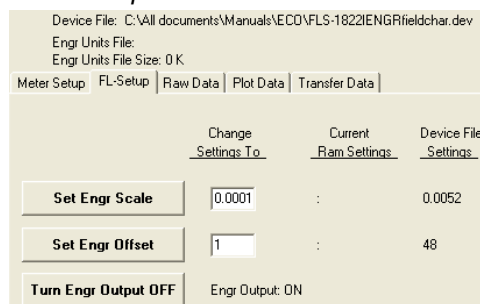

6. Enter the Scale Factor from the field characterization in the variable box under the Change Settings To column. This is the same value as in the edited device file.
7. Push **Set Engr Scale**.
8. Enter the Offset from the field characterization in the variable box under the Change Settings To column.
9. Push **Set Engr Offset**.
10. Push **Store to Flash**.

The sensor stores the field characterization values in its memory. The values show under the Current Ram Settings column in the host software.

3.5 Chlorophyll fluorometer-specific setup

ECO sensors that measure only chlorophyll have two device files. One is a standard device file. The other has an additional column that the host software uses to give chlorophyll output in $\mu\text{g/L}$.

1. Edit column 5 of the IENGR device file to show field characterization values.

```

ECO      FL-784
Created on:    07/17/11

:      ch1=ug/l
:      iengrunits =  $\mu\text{g/l}$  for CHL. ppb for PC, PE, CDOM, uranine.
:      column 5 = input scale factor and offset.

maxvoltage=   4.96
asv1=         6.3834
asv2=         12.7597
asv4=         25.5050

: Has internal CHL in meter output

COLUMNS=6
N/U=1
N/U=2
IENGR=3
N/U=4
Ch1=5   0.0077  81
N/U=6

```

2. Save the device file with a new name.
3. Refer to the steps in [Store field characterization values in sensor](#) on page 14 to save these values to the sensor in addition to the device file.
4. Select the *File* menu, then "Load Device File" in the host software.
5. Select the newly named device file.
6. Go to the *FL-Setup* tab.
7. Push **Turn Engr Output ON** to turn the $\mu\text{g/L}$ output on. Push **Start Data**. An additional column of data in units of $\mu\text{g/L}$ shows in the *Raw Data* tab.

Meter Setup	FL-Setup	Raw Data	Plot Data	Transf
02/06/12 10:17:40	28.22	695	3693	553
02/06/12 10:17:41	28.20	695	3691	553
02/06/12 10:17:42	28.23	695	3694	553

3.6 Device files

The host software uses a sensor-specific device file to show the data in the *Plot Data* tab and to calculate the data output in engineering units. Each device file has three necessary elements. The device file is not necessary to use the host software to set up and move data from a sensor.

1. The heading for the *Plot Data* tab.
2. The number of columns in the device file.
3. A description of the contents in each column.

Plot Data tab heading

The first line of the device file shows the model number and the serial number of the sensor. This information appears on the top of the *Plot Data* tab in the host software.

Number of columns

The column count shows how many columns of data the host software will process. The format is COLUMNS=x.

Description of columns

Each column of the sensor's data output comes from a description in the device file.

Number of columns=x

Date=x MM/DD/YY

Reference topics

Time=x HH:MM:SS
 N/U=x not used
 sc=scale factor
 off=offset
 IENGR=x
 mw=measurement wavelength of the sensor
 dw=display wavelength of the sensor

Example of a fluorometer device file	
chl, phycoerythrin, phycocyanin, uranine, rhodamine, or CDOM=x sc off Column 4 = scale factor (sc) offset (off).	<pre> ECO FLS-1822 Created on: 04/29/2011 : chl=ug/l. : iengrunits=ug/l for chl; ppb for PC, PE CDOM and uranine. : column 4=input the scale factor and offset in this column. : N/U=not used. maxvoltage=4.96 asv1=6.2606 asv2=12.5355 asv4=25.2860 COLUMNS=5 N/U=1 N/U=2 N/U=3 Chl=4 0.0052 48 N/U=5 </pre>
Example of a fluorometer device file with µg/L	
chl, phycoerythrin, phycocyanin, uranine, rhodamine, or CDOM=x sc off Column 5 = scale factor (sc) offset (off).	<pre> ECO FL-784 Created on: 07/17/12 : has internal CHL in output : iengr units = ug/l for CHL. ppb for PC, PE, CDOM, uranine. : column 5 = input the scale factor and offset values. maxvoltage=4.96 asv1= 6.3834 asv2= 12.7597 asv4= 25.5050 COLUMNS=6 N/U=1 N/U=2 IENGR=3 N/U=4 Chl=5 0.0052 48 N/U=6 </pre>
Example of a scattering sensor device file	
lambda (scattering wavelength) = x sc off mw dw Column 4 = scale factor (sc), offset (off), measurement wavelength (mw), and display wavelength (dw).	<pre> ECO BBS-974g Created on: 08/28/12 Columns=5 Date=1 Time=2 N/U=3 Lambda=4 7.916E-06 51 532 532 N/U=5 </pre>
Example of a turbidity sensor device file	
NTU = x sc off Column 4 = scale factor (sc) and offset (off).	<pre> ECO NTUSB-503 Created on: 09/07/2012 COLUMNS=5 N/U=1 N/U=2 N/U=3 NTU=4 0.0153 50 N/U=5 </pre>

3.7 Terminal program operation

Use Windows HyperTerminal® or other terminal program to operate sensors as an alternative to the host software.

Interface settings				
baud rate: 19200	stop bits: 1	data bits: 8	flow control: none	parity: none

3.7.1 Common terminal program operations

Command	Parameters	Description
!!!!	none	Stops data collected by the sensor. Allows the user to enter setup values. If the sensor is in a low-power state, turn the power supply off for one minute, then turn the power on and push the "!" key 5 or more times.
\$save	1–65535	The number of measurements that make up each row of collected data.
\$mnu	—	Prints the menu of available setup values to the host PC screen.
\$pkt	0–65535	Sets the number of rows of data that are collected between the specified time intervals.
\$run	—	Uses the current setup values to operate.
\$sto	—	Stores the desired setup values to the sensor's flash memory.

Single-parameter sensors—Fluorometer and NTU only

\$asv	1 2 4	Sets the sensor's analog scaling value. 1 = the analog data collected covers the bottom quarter of the sensor's data range. 2 = the analog data collected covers half of the sensor's data range. 4 = the analog data collected covers the sensor's entire data range.
-------	-------------	---

Fluorometer-only


\$scal	1 = ON 0 = OFF	Turns on the column with engineering units that show in µg/L. Turns off the column with engineering units that show in µg/L.
\$ugl	0–255	Sets the scale factor for data collected that shows in µg/L.
\$off	0–255	Sets the offset for data collected that shows in µg/L.

3.7.2 Terminal program operation for sensors with internal memory

Command	Parameters	Description
\$clk	24-hour time	Sets the time in the internal memory in the format hhmmss.
\$date	date	Sets the date in the internal memory in the format mmddyy.
\$emc	—	Clears the internal memory.
\$get	—	Reads data from the internal memory. Prints etx when it is complete.
\$int	24-hour time	Sets the time interval between sets of measurements the format hhmmss.
\$mvs	1 = ON; 0 = OFF	1 = the Bio-wiper is open. 0 = the Bio-wiper is closed.
\$rec	1 = ON 0 = OFF	1 = Turns on the sensor's internal memory. 0 = Turns off the sensor's internal memory.
\$rls	—	Loads the settings from the flash memory.
\$set	0–65535	Sets the number of rows of data that are output between low-power states.

3.7.3 ECOView and terminal program differences

ECOView uses different but equivalent words for data collection options.

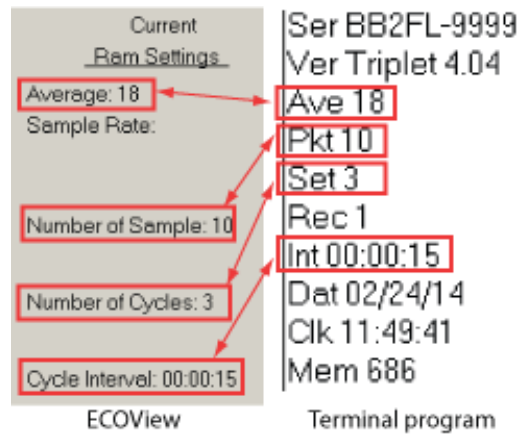

Section 4 Optional equipment

4.1 Test cable

Use a test cable to set up and test the sensor before deployment.



One analog output connector



Two analog output connectors



No analog output

1 six-contact connector	3 db-9 serial port connector
2 9-volt battery connector	4 RCA connector

1. Connect the six-contact connector into the sensor.
2. Connect the 9-volt connector to a 9-volt battery. As an alternative, it can be connected to a regulated power supply.
3. Connect the db-9 connector to the host PC. Use a USB-to-RS232 adapter cable if necessary.
4. Use a digital multimeter (DMM) to see analog output from the sensor. The inside of the RCA is signal (red DMM probe) and the outside is ground (black DMM probe).

4.2 Internal batteries

Six 9-volt lithium batteries supply power to sensors that have batteries. The sensors can also use alkaline or lithium manganese dioxide (LiMnO₂) batteries. Alkaline batteries supply about 1000 mA-hours of power. LiMnO₂ batteries supply more than 2000 mA-hours of power.

Note: *The nominal water temperature, sequence timing, sample periods, and other variables will change the use time of the sensor's internal batteries.*

4.2.1 Replace batteries

⚠ WARNING

The pressure housing of the *ECO* sensor must be opened to replace the batteries. If done improperly, it could result in personal injury or death due to abnormal internal pressure as a result of flooding. It may not be possible to repair flooded sensors.

The manufacturer disclaims all product liability from the use or service of this equipment. The manufacturer cannot control the use of these sensors or choose qualified personnel to operate them, and therefore cannot take steps to comply with laws pertaining to product liability, including laws that impose a duty to warn the user of any dangers involved with the operation and maintenance of the sensors. Acceptance of this equipment by the customer shall be conclusively deemed to include a covenant by the customer to defend and hold the manufacturer harmless from all product liability claims arising from the use and servicing of this equipment. Flooded instruments will be serviced at the discretion of the manufacturer.

Flooded sensors are usually not serviceable. There is a small chance that the manufacturer can get the data that is stored in the sensor. Contact service@wetlabs.com for more information about flooded sensors.

1. Clean any debris from the end flange.
2. Dry the sensor thoroughly.

3. Remove the dummy plugs if necessary.

⚠ CAUTION

The sensor may be under pressure. Do not point it toward any part of the body when removing the vent plug or end flange.

⚠ CAUTION

Replace the batteries in a clean and dry environment. Gases in the sensor can expand and push the pressure relief plug open. This will flood the sensor. Do not replace batteries in a cold environment and then deploy the sensor in a hot environment.

4. Point the connector end flange down, away from the face.
 - a. Pull the vent plug loose.
 - b. If the sensor has an external thermistor, pull it loose.
5. Dry the vent plug (and the thermistor, if applicable).
6. Use needle-nose pliers to remove the monofilament from the end flange.

Figure 3 Pull monofilament from end flange




7. Remove the end flange from the pressure housing.
The jacking screws supplied as spare parts can be used to push the end flange off of the pressure housing, and then be removed.
8. Gently disconnect each Molex[®] connector.
9. Remove the screw that holds the vent plug in the end flange.
10. Dry the end flange and the pressure housing seal areas.
11. Examine the O-rings on the vent plug and thermistor (if applicable).
Remove any O-ring that is damaged.
12. Apply a light coat of vacuum grease on a new O10 O-ring and put it on the vent plug or thermistor.
13. Put the vent plug into the top of the end flange.
14. If applicable, insert the thermistor into the end flange again.
15. Put the screw for the vent plug in the inside of the end flange.
The screw holds the vent plug in the end flange.
16. Carefully pull the white plastic loop to remove the battery pack from the pressure housing.
17. Remove the black plastic protectors from the ends of the long screws securing the batteries.
18. Loosen, but do not remove, the retaining screws using a 1/4-inch slotted driver.
Do not remove both screws. Removing both of the retaining screws will result in a pile of parts and will complicate the replacement process.

Figure 4 Both retaining screws removed




19. Disconnect each of the six batteries.
20. Tilt the battery board enough to connect the first battery in the contacts that are perpendicular to the other two.

Figure 5 Connect the first battery

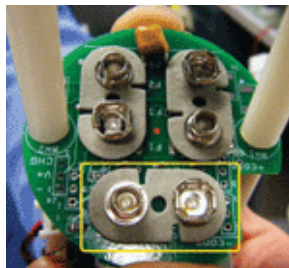

21. Pivot the boards in the opposite direction to connect the other two batteries.
22. Connect the second set of the batteries.
23. Hold the assembly plates on top and bottom and tighten the screws.
The bottoms of the batteries may splay out. Make sure that the batteries do not stick out beyond the circuit boards. If they do, the battery pack will scratch the sealing surface when it is put into the pressure housing again.
24. Install the bottom neoprene pad and the black plastic protective covers on the ends of the screws.
25. Remove and examine the 224 O-ring on the pressure housing for damage.
26. Use a new O-ring if necessary.
27. Apply a light layer of grease such as Dow Corning® High Vacuum Grease on the O-ring.
28. Put the battery pack into the pressure housing.
29. Attach the Molex® connectors.
Note that the dowel pin may be in the end flange or in the pressure housing.

Figure 6 Inside of end flange

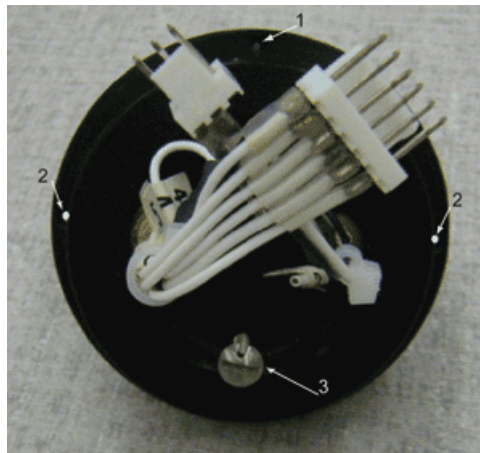

1 dowel hole

2 through-holes for jacking screws

3 vent plug screw

30. Put the end flange onto the housing so that the wires are out of the way.
31. Align the dowel pin with the dowel hole in the end flange, not the jack screw holes, which go through the end flange.
32. Make sure that the wires are not pinched between the end flange and the pressure housing.
33. Make sure that the end flange is attached to the pressure housing.
34. Install the monofilament in the end flange.

4.3 Bio-wiper and copper faceplate

The Bio-wiper and copper faceplate help decrease bio-fouling during extended deployments. The Bio-wiper is controlled manually using a host controller, or it can work autonomously as pre-programmed before a deployment. The temperature and the depth at which the sensor is deployed will change the Bio-wiper cycle time.

4.3.1 Clean Bio-wiper and faceplate


Sensors that have copper faceplates and copper Bio-wipers need to be cleaned at regular intervals to keep the copper effective as an anti-fouling aid. Remove the Bio-wiper and faceplate from the sensor to clean them.

NOTICE

Do not twist the Bio-wiper or the Bio-wiper shaft by hand. Push **Open Shutter** and **Close Shutter** in the host software to turn the Bio-wiper. Turning by hand will damage the Bio-wiper motor.

1. Disconnect the sensor from any power source.
2. Use the manufacturer-supplied 3/32-inch hex key to loosen the screw that holds the Bio-wiper to the sensor.
If necessary, insert the manufacturer-supplied jacking screw into the expander hole to loosen the Bio-wiper. *PAR sensors do not have an expander hole and do not ship with jacking screws.*


Figure 7 Bio-wiper screw holes

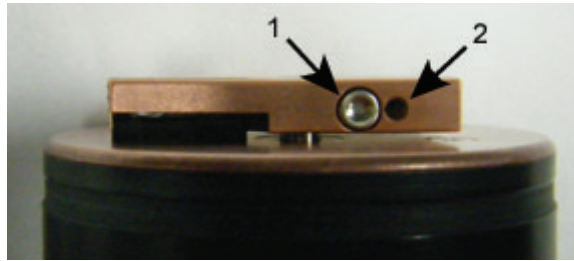

1 the screw in the clamping hole	2 the expander hole
----------------------------------	---------------------

3. Lift the Bio-wiper off of the faceplate.
4. If used, remove the jacking screw.
5. Use a small Phillips screwdriver to remove the screws that attach the faceplate to the optics face.
6. Keep the screws.
7. Clean the Bio-wiper and faceplate with soapy water.
8. Rinse and dry thoroughly.
9. Use a Scotch-Brite® or similar pad to buff the Bio-wiper and faceplate until they are shiny.
10. Use an isopropyl alcohol-saturated cotton swab to clean the Bio-wiper shaft and the shaft hole.
11. Make sure the Bio-wiper and faceplate dry completely.
12. Attach the faceplate again.
13. Make sure the screw used to attach the Bio-wiper is in good condition.
The 3/32-inch hex key must fit into the screw's socket.
14. If the screw is damaged, use a new screw from the spare parts kit that the manufacturer supplies with the sensor.
These screws are 4-40 x 3/8-inch 316 stainless steel and are treated with anti-seize.
15. Insert a jacking screw into the expander hole of the Bio-wiper. Turn it slowly until the Bio-wiper goes onto the shaft easily.
16. Put the clean Bio-wiper over the shaft but do not twist the Bio-wiper.
17. Make sure to remove the jacking screw before tightening the 3/32-inch clamping screw.
18. Use the commands in the host software to turn the Bio-wiper into the closed position.
The Bio-wiper covers the optics.
19. Set the clearance between the Bio-wiper and the faceplate to approximately 0.8 mm.
Note: If the Bio-wiper is has too much bend because it is too close to the faceplate, the motor will use too much current. If it is not close enough to the faceplate, the Bio-wiper will not keep the optics clean.

Figure 8 Bio-wiper clearance and wiper bend

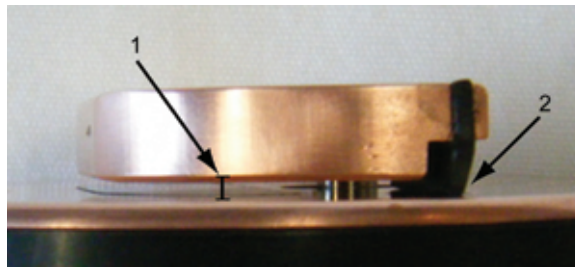

1 Bio-wiper clearance from faceplate


2 Bio-wiper flex

20. Use the 3/32-inch hex key to tighten the screw. Do not over-tighten the screw.
21. Turn on the sensor to make sure that the Bio-wiper rotates correctly.
The Bio-wiper must rotate 180 degrees and uncover the optics before they come on.
The Bio-wiper will then rotate 180 degrees to cover the optics after they come on.
22. Make sure that the Bio-wiper operates correctly.

4.4 External thermistor

The calibration coefficient value for the thermistor is on the sensor's characterization page that ships with the sensor. The thermistor gives temperature output in counts. Do one of the methods below to change counts to engineering units.

1. Use the drop-down menu in the host software to see the thermistor output in °C.

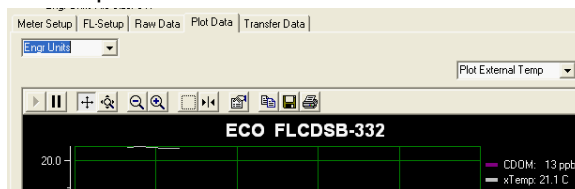


2. Use MATLAB, MS Excel or other software to solve for:
Temperature, °C = (Output × Slope) + Intercept

4.5 Pressure sensor

The calibration coefficient value for the pressure sensor is on the calibration page that ships with the sensor. The pressure sensor gives pressure output in counts. Do one of the methods below to change counts to engineering units.

1. Use the drop-down menu in the host software to see the pressure sensor output in dbar.

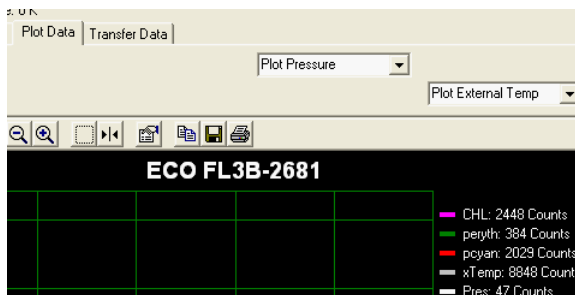

2. Use MATLAB, MS Excel or other software to solve for:
relative pressure, dbar = (output × slope) + intercept
3. Calculate the absolute pressure:
absolute pressure, dbar = relative pressure, dbar - relative pressure at atmospheric-water interface, dbar


Pressure sensors need to be set at zero for each deployment. Do not deploy the pressure sensor deeper than the depth rating on the calibration page.

4.5.1 Maintain pressure sensor

The plastic fitting filled with silicone oil is a buffer between the pressure transducer diaphragm and seawater. Add silicone oil to the reservoir on top of the transducer at regular intervals.

1. Make sure the top of the sensor is clean.
2. Use a 9/16-inch wrench to hold the white nylon Swagelok® fitting.
3. Use a 7/16-inch wrench to loosen the cap on top of the fitting.
4. Remove the cap.
5. Use a wire or toothpick to clean the hole in the cap.
Do not blow compressed air into the fitting. It will make a mess.
6. Add Dow Corning® 200 Silicone Oil to the reservoir until the oil is visible.

Figure 9 Pressure sensor cap

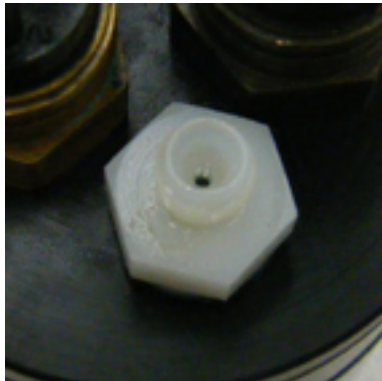

7. Replace the cap. Make sure not to tighten the cap too much.
8. Clean any excess oil from the end flange of the sensor.

Section 5 General information

Revised editions of this user manual are on the manufacturer's website.

5.1 Warranty

This sensor is warranted against defects in materials and workmanship for one year from the date of purchase. The warranty is void if the manufacturer finds the sensor was abused or neglected beyond the normal wear and tear of deployment.

5.2 Service and repair

The manufacturer recommends that sensors be sent back to the factory annually to be cleaned, calibrated, and for standard maintenance. Do the steps below to send a sensor back to the manufacturer.

1. Contact the manufacturer for a Return Merchandise Authorization (RMA).
Note: The manufacturer is not responsible for damage to the sensor during return shipment.
2. Remove all anti-fouling treatment from the sensor before sending it back to the manufacturer.
Note: The manufacturer will not accept sensors that have been treated with anti-fouling compounds for service or repair. This includes tri-butyl tin, marine anti-fouling paint, abrasive coatings, etc.
3. Use the sensor's original ruggedized shipping case to send it back to the manufacturer.
4. Write the RMA number on the outside of the shipping case and on the packing list.
5. Use 3rd-day air to ship the sensor back to the manufacturer. Do not use ground shipping.
6. The manufacturer will supply all replacement parts and labor and pay to send the sensor back to the user via 3rd-day air shipping.

5.3 Support email addresses

Refer to the website for FAQs and technical notes, or contact the manufacturer for support at:

support@wetlabs.com




5.4 Waste electrical and electronic equipment



Electrical equipment that is marked with this symbol may not be disposed of in European public disposal systems. In conformity with EU Directive 2002/98/EC, European electrical equipment users must return old or end-of-life equipment to the manufacturer for disposal at no charge to the user. To recycle, please contact the manufacturer for instructions on how to return end-of-life equipment, manufacturer-supplied electrical accessories, and auxiliary items for proper disposal.

Appendix A *ECO* mounting bracket

The details of the mounting bracket for the *ECO* sensors are shown below. The mounting bracket is not installed on sensors with internal batteries or those rated to 6000 m depth.

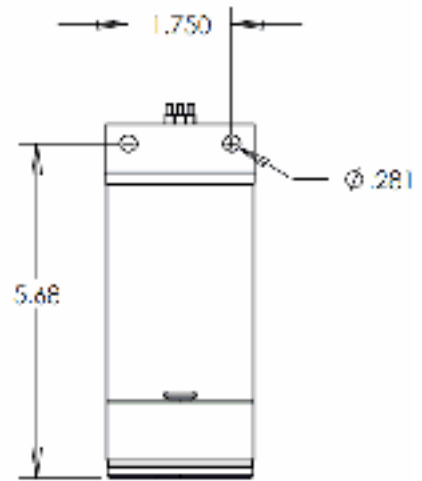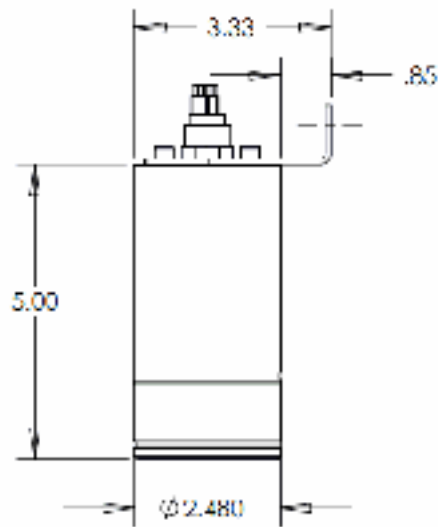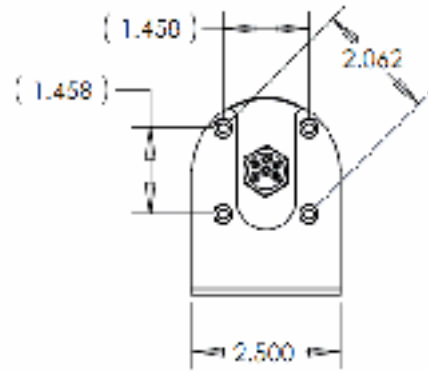

WET Labs, Inc.
P.O. Box 518
Philomath, OR 97370 U.S.A.
Tel. (541) 929-5650
Fax (541) 929-5277
www.wetlabs.com

